


# RING LIGHT WITH HOLDER for Rapidtookit tests

A low cost alternative to fiber optic cold lighting or LED ring lights, developed specifically for the Rapidtookit.

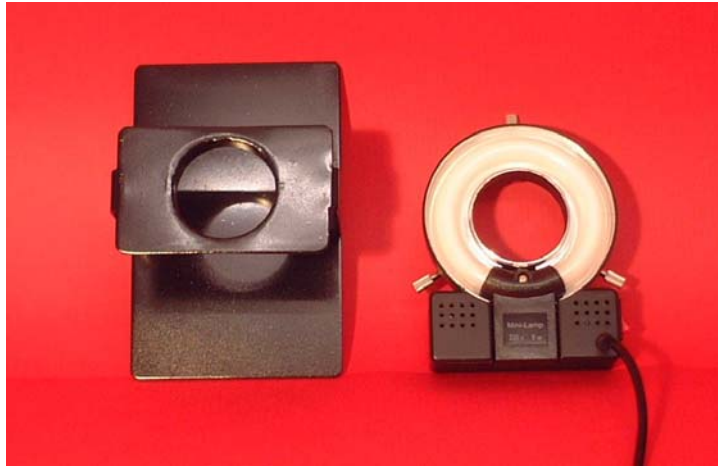

## RING LIGHT

- Fluorescent circular strip light (220V–8W, also available in 110V) of 8.5 cm diameter, mounted into a circular metal frame provided with an on-off switch.
- The ring light fits directly on the bottom of the stereohead of some types of microscopes (but not on all types).
- A second strip light is provided as spare part.

## HOLDER

- A unique holder in stainless steel which can be placed on the bottom stage of any type of stereomicroscope.
- The holder consists of a flat top plate (6 x 10 cm) with a central tube of 4.5 cm diameter, positioned 7 cm above the (16 x 10 cm) bottom plate.
- The ring light can be fixed on or removed very easily from the central tube of the holder with the 3 screws located on the outside of the metal frame of the ring light.